

About head lice

Every case of lice begins with one louse or a few lice located on the head, laying eggs that hatch within 7-10 days and continue the 3-stage cycle of reproduction and infestation:

Lice eggs: Usually laid close to the scalp, on the back of the head and behind the ears. Once the louse has hatched from the egg (within 7-10 days), the empty case remains stuck to the hair.

Nymphs : Young lice that have hatched from the eggs. At this stage, they are not capable of reproducing.

Mature lice: Lice become mature 10 days after hatching and begin to reproduce after reaching full size, at which stage egg laying begins. The size of head lice varies according to their age, ranging from 0.5-3 mm in length.

Head lice can be transmitted only through direct head-to-head contact. They hold on to the hair with claws, but after treatment they lose their strength and fall off.

1 Before you start

Please read this user manual in full before using your Novokid Lice Treatment Kit. Be sure to read all warnings before use, and take all the recommended safety precautions.

1.1. Device description

Novokid is an innovative device for the effective treatment head lice, super lice and eggs.

Novokid uses a natural, plant-based solution that is non-toxic, free of pesticides, silicon or alcohol, efficient and easy to use.

1.2. Package contents

Your Novokid Lice Treatment kit contains: 1 Novokid vaporizer device
4 Novokid treatment capsules 1 Novokid cap
2 Connector pins
1 Silicone tube
1 Lice comb
1 Screwdriver for battery replacement
AA Batteries – 4 pack
1 User guide & warranty certificate

2 Intended use and safety

Novokid has been scientifically proven to be suitable for use on all types and lengths of hair, for use on children aged 4 and up, under adult supervision. Do not stand behind the child's head while the device is on. Novokid Lice Treatment is for external use only.

3 Warning

Do not use in the case of infection, wound, skin irritation, or skin problems. If a skin rash or any other sign of allergic reaction develops during or after use, discontinue treatment. Wash the skin with water immediately. If irritation persists, seek medical advice.

4 Possible side effects

The treatment may cause itching or a slight irritation on the scalp.

5 Before treatment

The treatment may cause itching or a slight irritation on the scalp.

6 Get started!

The Novokid Lice Treatment takes 10 minutes. Follow these simple steps:
Place the cap on the child's head so that the "Back" sticker is on the back side of the child's head and that the white band lies flat across the forehead.
Make sure that the hair and ears are fully covered by the cap.
Do not cover the eyes with the cap.
Be sure to gently squeeze the air out of the cap. This is critical for the treatment's success.
Take a single capsule and click it into connector placed on the cap. Connect the other end of the capsule to the end of the tube - when it is attached securely you will hear a click.

You are now ready to start the treatment!

The treatment may cause itching or a slight irritation on the scalp.

Notes

The Novokid Lice Treatment takes 10 minutes. Follow these simple steps:
Place the cap on the child's head so that the "Back" sticker is on the back side of the child's head and that the white band lies flat across the forehead.
Make sure that the hair and ears are fully covered by the cap.
Do not cover the eyes with the cap.

Be sure to gently squeeze the air out of the cap. This is critical for the treatment's success. Take a single capsule and click it into connector placed on the cap. Connect the other end of the capsule to the end of the tube - when it is attached securely you will hear a click.

7 Storage and maintenance

Please read this user manual in full before using your Novokid Lice Treatment Kit. Be sure to read all warnings before use, and take all the recommended safety precautions.

8 Consumables

We recommend that:

Keep your Novokid Lice Treatment out of the reach of children.
Keep all parts in the box, and store at room temperature for future use. you Do not wrap the tube around the device.
Do not use a capsule after the expiration date on the packet. Replace batteries once the red light appears on the device.
Use soap and water to clean the cap, and wipe it using a damp cloth
Do not allow the base unit to become wet.

9 Troubleshooting

If the cap does not inflate once the device is switched on, tighten the capsule connection at the back of the head and check that the relief valve is in place. If the LED light does not flash, check that the batteries are positioned correctly, or replace them. For all other problems, please contact Novokid Customer Service.

10 Customer service

See our website for more information or for customer support: www.Techcareltd.com. Customer service

11 Warranty

Your Novokid Lice Treatment kit is covered by a 3 Year warranty. Please refer to your Novokid warranty card for further details.

12 Disposal

Please dispose of packaging at a recycling point that sorts materials by type. Dispose of used batteries separately from household waste and in accordance with your country's waste regulations.



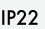








13 Specifications

Your Novokid Lice Treatment kit is covered by a 3 Year warranty. Please refer to your Novokid warranty card for further details.

Model	TC-001
Power requirements	AA batteries x 6V = 500MA 4
Operating temperture & humidity	:C 15-30 Humidity 70%>
Pump pressure	LPM 4
Storage instructions	Room temperature (See paragraph 7)

14 Signs and symbols

Your Novokid Lice Treatment kit is covered by a 3 Year warranty. Please refer to your Novokid warranty card for further details.

Products marked with this symbol meet the requirements of the CE directive.	
For indoor use only	
The device is protected against the insertion of objects bigger than 2.5 mm diameter, such as fingers, to protect against electric shock.	
Do not use in a wet environment.	
Do not dispose with household waste. Follow your country's rules for the separate collection of electrical and electronic products and batteries. Correct disposal helps prevent negative consequences for the environment and human health.	
Operating temperature range 15-30 °C.	
The user manual contains important warnings and instructions for use.--	
Date of manufacture	
Authorized representative in the European	
Use by date	
NOVOMIC LTD. A subsidiary of TechCare corp LTD. 23 Hamelacha St. Rosh Haain, Israel	
MEDES LTD.5 Beaumont Gate, Shanley Hill. Radlett Herts WD7 7AR, United Kingdom. Tel/Fax : +44(0)1923859810	